
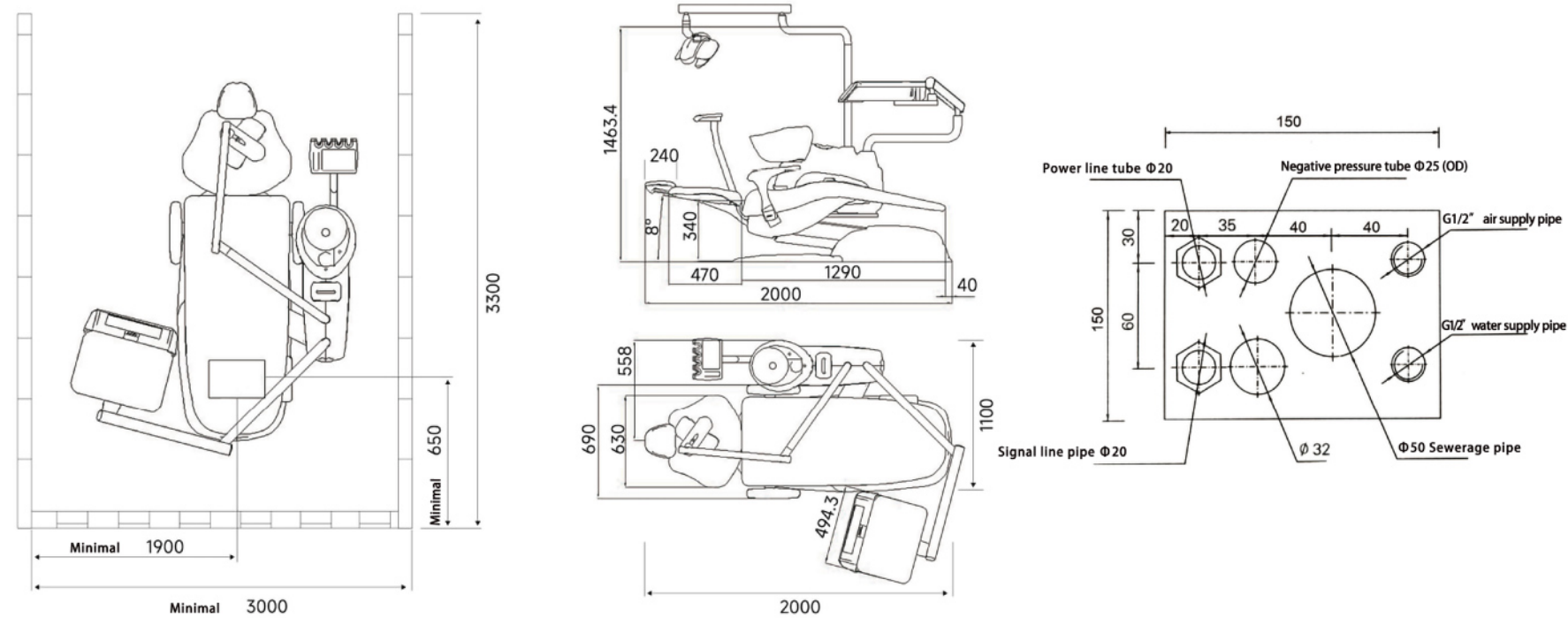


1. water supply pipe: white PPR hot melt pipe, reserved G1/2 "copper internal tooth interface (recommended to leave G1/2 "internal and external teeth through the valve);
2. Air supply pipe: seamless copper pipe or stainless steel pipe, reserved G1/2" copper inner tooth interface (recommended to leave G1/2" inner and outer teeth through the valve);
3. Drainage pipe: white PVC drainage pipe, main pipe Φ75mm, branch pipe Φ50mm.
4. Transmission line: with PVC Φ32mm cable sets HD SDI-75-5 line (BNC connector), more than 50 meters of transmission distance is recommended to use fiber optic HDMI cable;
5. Power line: white PVC wiring pipe Φ20mm, tube sets 3 * 1.5m line, fire, zero, ground line from the air switch box, reserved with a switch outlet;
6. Negative pressure signal line: white PVC pipe, Φ20mm tube set 2 * 1m line.
7. Negative pressure pipe: blue PVC pipe main Φ50mm, branch Φ32mm, access to dental chair pipe Φ25mm.

Note: Drainage pipe and negative pressure pump receiver can not be used at right angles to the joints, the application of 45-degree elbow. Purifier water supply and air compressor into the air should be reserved switch.

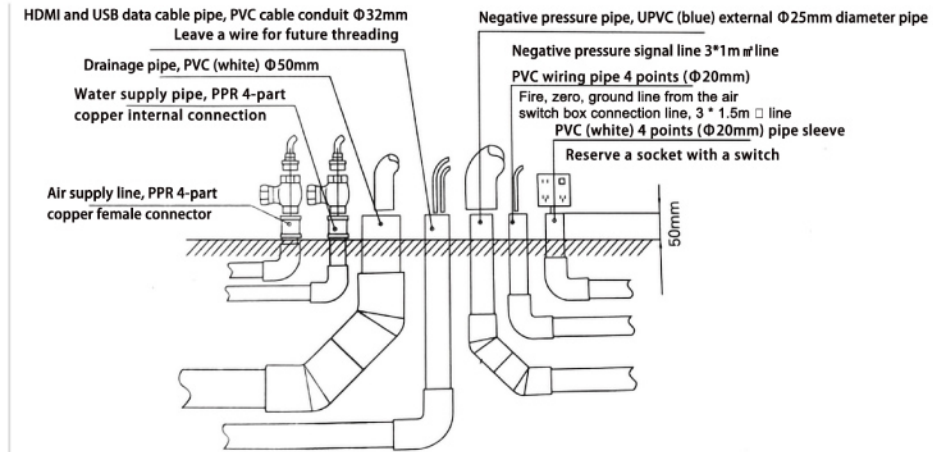

website: www.dentalstore-eg.com

Contact number: 01070873434

Dental Unit

Mirage Q9 Dental Unit



Digital code touch control

· Digital display and touch control, perfectly adhering to ergonomic design, enable doctors to precisely control and have a clear understanding during the operation process.



Dual-joint multi-angle adjustable assistant stent

· The flexible rotation design allows for free adjustment of the Angle according to the patient's position




Rotatable ceramic spittoon

· Ceramics, automatic constant temperature heating system



Flexible foot pads

· 130mm forward extension

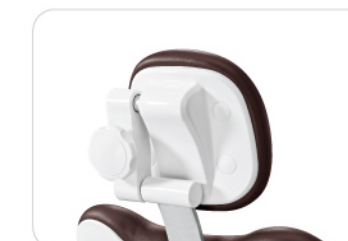

Multi-dimensional adjustable headrest

· Fine-tuning the Angle and sliding forward and backward precisely support the neck curve

Professional Support The New Standard Of Comfort

Detachable neck pillow lumbar support assembly

· Modular design, easy disassembly, and convenient for cleaning and disinfection



Mirage Q8 Dental Unit



Aromatherapy
Soothes anxiety

As Steady As A Rock Silently Guarding



Metal backrest

·Sturdy and durable, it is not prone to deformation even after long-term use




Fragrant

·The fragrance can release a relaxed aroma



Multifunctional foot switch

·Lift, water and gas, mouthwash, rush paddy

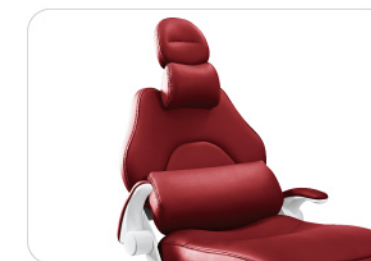


Detachable neck pillow lumbar support assembly

·Modular design, easy disassembly, and convenient for cleaning and disinfection



UVC water disinfection device

·The molecular structure of DNA or RNA that directly kills harmful microorganisms in the cells of the organism.

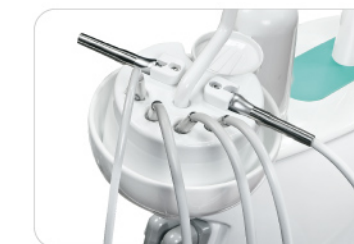

Automatic disinfection system

·Clear biofilms in the pipeline, eliminate bacteria and viruses


Mirage Q7 Dental Unit



Aromatherapy
Soothes anxiety



GoMed
Medical Solutions

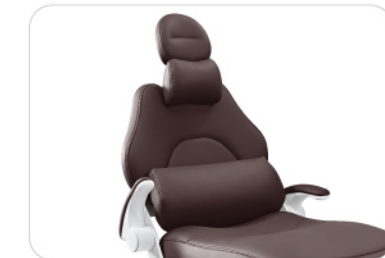


Detachable neck pillow lumbar support assembly
· Modular design, easy disassembly, and convenient for cleaning and disinfection



UVC water disinfection device
· The molecular structure of DNA or RNA that directly kills harmful microorganisms in the cells of the organism.

Meticulous Attention To Detail
Quality Speaks For Itself

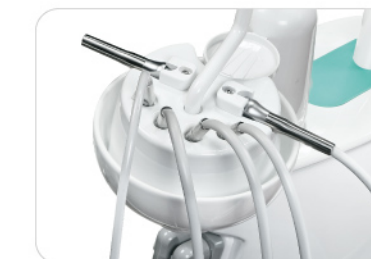

One-touch automatic disinfection system
· Clear biofilms, bacteria and viruses in the pipeline



Fragrant
The fragrance can release a relaxed aroma



Multifunctional foot switch
· Lift, water and gas, mouthwash, rush paddy

Mirage Q6 Dental Unit



GoMed
Medical Solutions

LED Oral lamp(Six -hole yellow /white light)
· Inductive switch control, stepless brightness adjustment

As Steady As A Rock
Silently Guarding



Rotating box design
· Doctors can operate from multiple angles to reduce the number of adjustments made by patients



Disassembled saliva filtering system
· Prevent pollutants from entering the pipeline system

Rotatable ceramic spittoon
· Ceramics, automatic constant temperature heating system



Optional Customization



Skin sensation color



PU cushion color



Product standard Configuration

CONFIGURATION	MODEL NO.			
CONFIGURATION NAME	Mirage Q9	Mirage Q8 One-click disinfection	Mirage Q7 One-click disinfection	Mirage Q6 Updated Version
24v DC Motor	Timotion	Timotion	Timotion	Domestic
3-way Syring 2 Pic(Hot/cold)	●	●	●	●
Main Control/ Uxiliary Control	●	●	●	●
Flushing/water Quantify Supply System	●	●	●	●
Water Supply Heating System	●	●	●	●
Aluminum Alloy Strong/weak Salivary Ejector	●	●	●	●
Working Disc With Air Brake	●	●	●	●
Protective Film Of Instrument Tray	SilicaGel	○	○	SilicaGel
Instrument Tray With Stainless Steel Tray	○	●	●	○
X-ray Viwer	●	●	●	●
Purified Water/tap Water Conversion System	●	●	●	●
90 ° rotating ceramic spittoon	●	●	●	●
Multi-function Foot Switch(Adjustable)	●	Luxury	Luxury	●
The Cushion	Skin sensation	Skin sensation	Skin sensation	Skin sensation
Oral Lamp	6 holes	8 holes	8 holes	6 holes
Position With Memory	9 Postiton	9 Postiton	9 Postiton	○
Geramic Trolley	○	○	○	○
Water/air/electricity Main Switch	●	●	●	○
Solenoid Valve	●	●	●	●
Armrest	Double	Double	Double	Double
Store Water Bottles	●	Double (sterilized/purified)	Double (sterilized/purified)	●
Doctor Chairs	K-03	K-03	K-03	K-04
Cast Aluminum Base	14mm Steel plate	14mm Steel plate	14mm Steel plate	14mm Steel plate
Air/water Pipe	Imported	Imported disinfection	Imported disinfection	Imported
Reserve Interface For Scaler	●	●	●	●
Tissue Box	●	●	●	○
Handpiece water heating	○	●	●	○
Options				
Endoscope Camera	○	○	○	○
Holder Of Endoscope Camera	○	○	○	○
EGG Monitor	○	○	○	○
Luxury Shadowless Operating Lamp	○	○	○	○
Up-mounted Instrument Tray	○	○	○	○

MARK: ● Configure ○ Optional